



ivo Lab

Ivostar® | Gnathostar®

Simply good denture teeth



Well-proven
and economical
– for satisfied
patients



Come see us
at the 2013 DENTAL
EXPO – INTERACTIVE
17 & 18 May 2013
www.dentalexpo.co.nz

Inexpensive dentures

- Wide range of indications
- Tried-and-tested material





Ivostar | Gnathostar

Simply good denture teeth

// Talk to us about consignment or trade-in options //



Need a wider range of moulds?

- Extensive choice of moulds available to suit each patient
- Mould guides available

Denture teeth debonding?

- Special acrylic in neck area ensures strong bond to denture base
- Proven PMMA construction

Looking for a lower priced toothline, that does not compromise aesthetics or quality?

- Excellent aesthetics
- 3 layer anterior / 2 layer posterior
- Full shade range – 16 Vita A-D and 4 BL Bleach
- Used for over 10 years by denture wearers

Need teeth now?

- We carry the stock so YOU can have your teeth on hand when you need them – consignment options available





Added value for your laboratory at a reduced cost!

Introducing... **Harmony[®] 2**

Crown & Bridge Alloy

NEW

Our new, most economically efficient yellow, palladium-silver-based crown and bridge alloy

- Onlays, $\frac{3}{4}$ crowns, crowns, telescope and conus crowns, root canal posts, short- and long-span bridges
- Tested and certified biocompatibility
- Excellent corrosion resistance – crowns will not tarnish
- Economical due to a very low density – 9.8
- Simple processing and polishing properties and a brilliant high gloss is easy to achieve
- Suitable for bonding systems and veneering composites, such as SR Link and SR Nexco
- Au 2.0%, Pd 32.7%, Ag 34.3%, In 28.0, Zn 2.6%
- Yellow colour, Type 3


For your added security, precision & peace-of-mind, your alloy will shortly be delivered in new 5g sealed packets

IPS[®] Natural Die Material

Die shade selection

In order to facilitate the reproduction of the desired tooth shade, the shade of the preparation is determined with the help of the IPS Natural Die Material shade guide. This enables the technician to fabricate a model die similar to the preparation of the patient, on the basis of which the correct shade and brightness values of the all-ceramic restorations may be selected.

- ND1 – ND9 shades available

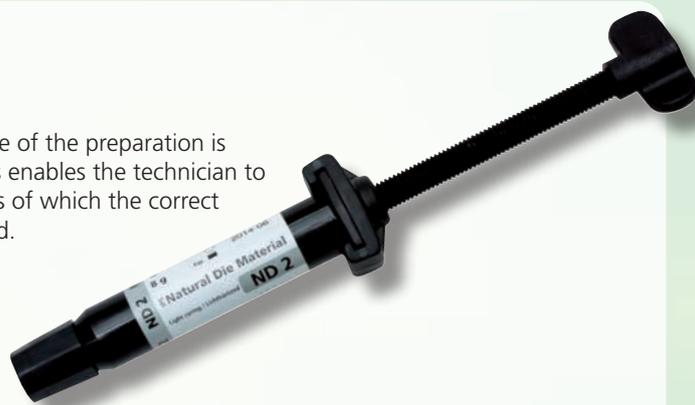

Example of the preparation shade effect

Crown made of IPS e.max Press HT B1 on different preparation shades.

IPS PressVEST Speed

For accurate press results with excellent crack resistance and the least amount of reaction layer. IPS PressVEST Speed is a phosphate-bonded investment material for the speed heating technique.



IPS e.max[®] Press Invex Liquid

The IPS e.max Press Invex Liquid is used to soften the surface reaction layer formed on pressed IPS e.max Press and ZirPress frameworks during the pressing procedure.





An ingot for EVERY indication...

IPS e.max® Press

Inlays, Onlays & Veneers

HT – THE CHAMELEON



Slightly discoloured preps

MO – THE BASE FOUNDATION



Occlusal veneers or no prep veneers

Impulse – THE ENAMEL



Full contour crowns

LT – THE ALL-ROUNDER



Severely discoloured preps, metal posts or implant abutments

HO – BLOCKER




Programat® EP3000/G2



Programat® EP5000/G2



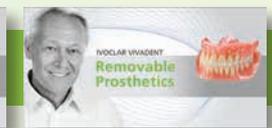

Stratos 100

The Average-Value Articulator



The Stratos 100 is a simple average-value articulator with very convenient handling. Ivoclar Vivadent teeth can be set-up easily and quickly in conjunction with the setting-up template. **Indication: Removable denture prosthetics, BPS System**

SIMPLE
PRACTICAL
UNCOMPLICATED

- Centric locking catch 
- 120° opening stop 
- Protrusion screw 
- 45° incline support holder 

Stratos 200

The Semi-Adjustable Articulator



The Stratos 200 is a semi-adjustable articulator with various setting options due to exchangeable joint inserts. It has been especially designed for use in fixed and removable denture prosthetics. **Indication: Removable denture prosthetics, BPS System, fixed crowns and bridges.**

VERSATILE
PRACTICAL
UNIVERSAL

- Exchangeable joint inserts 
- Centric locking catch 
- Protrusion screw 
- Side-shift setting 

Stratos 300

The Individual Articulator



The Stratos 300 is a state-of-the-art precision articulator with individual setting options. With its wealth of possible variations, it represents the optimum basis for exact functional simulations for demanding, high-quality dentures. **Indication: Removable denture prosthetics, BPS System, fixed crowns and bridges, implant-supported dentures, functional diagnostics**

HIGHLY ADJUSTABLE
VERSATILE
FLEXIBLE

- Bennett angle setting 
- Centric locking catch 
- Protrusive shift 
- 120° opening stop 
- Side-shift setting 



Telio Lab

Unlimited flexibility in the dental laboratory

Telio Lab is a powder/liquid-based resin for the fabrication of long-term temporary restorations using the cold technique.

Designing the shape with a wax-up and silicone key.



The material is polymerised under pressure in a water bath for 15 minutes at 2-6 bar and a temperature of 40-50 °C.

Telio Lab LC

Optimised application with supplementary materials

The Telio Lab LC components consist of light-curing paste-type incisal materials and flowable stains. They are designed for applying customisations and add-ons to Telio restorations.



Disposable Plaster Trap Bags

- Easy cleaning, less mess
- Installs into Plaster Trap neatly and lifts out to dispose

NEW



Plaster Traps

- Easy to install and clean plastic basins

VD/G (540 x 380 x 250mm)
VD/P (440 x 380 x 250mm)



Poliresin®

An odourless polishing agent specially developed for polishing acrylic dentures.

- Poliresin is used in the same way as pumice, but smoothes even more effectively and has a 30% longer life span
- The acrylic is gently polished due to the unique crystalline structure of Poliresin
- The time required for the subsequent high-lustre polish is greatly reduced
- EH 100323



GEO Crowax

- All-purpose wax – available in Green-transparent, Blue-transparent, Red-transparent, Beige-opaque, Grey-opaque
- R474



GEO Classic

- R497

GEO Avantgarde

- R495 - Universal
- R492 Occlusal



Giroform Die Link

Varnishes, Thinner & Hardener

A series of air-drying die varnishes which contain innovative chemical primer components as well as nano-particles which ensure excellent bonding to dental plasters or stones and provide smooth and scratch resistant surfaces. The superthixotropic consistency guarantees an even, homogeneous layer thickness even on sharp corners.

- Varnishes... Blue (20µm), Grey (20µm), Gold (15µm), Silver (13µm), Red (8µm), plus Hardener and Thinner
- 15ml bottle refills



* Special Price & Free Thinner ONLY while stocks last. 1 Free Thinner per invoice.

Hinrisol

Alcohol-Based Surface Wetting Agent

For reducing surface tension, neutralising and cleaning silicone, wax and resin surfaces without leaving a film. Hinrisol is compatible with all silicone-based impression and duplicating materials.



Hinrisid

Surfactant Debubbliser Also suitable for polyether!

Hinrisid is a universal wetting agent for fabricating gypsum models without bubbles in silicone and polyether impressions.



Hinridur® S

Synthetic stone Type 3

A creamy synthetic hard stone especially suitable for prosthetics and orthodontics. The super-white colour is specially recommended for orthodontics work, and the yellow for working models.

- Yellow & Super-White



All Orders over \$400
FREIGHT FREE*

*Exclusions Apply
Excludes GST

Courses & Events - 2013



Dates	Course/Event	Lecturer	Location
11-13 April	Implant Over-dentures + BPS Certification	Arnie Jabour	ICDE Auckland
16 - 17 May	NZACD Symposium	Various	Auckland
27 / 28 / 29 May	Aesthetic Challenge - Metal & Metal-Free	Dieter Grubel	Wellington, Christchurch, Auckland ICDE
30 May	The Secret of Aesthetics - Metal Ceramics	Dieter Grubel	Auckland ICDE
13 - 14 June	Advanced Techniques for Alloy & Cermaics	Jeff Smith	Auckland ICDE
27 - 28 June	Challenging Class II, Class III and Lingualised Occlusion	Steve McGlynn	Auckland ICDE

For detailed course information, contact Caroline Jones on (09) 914 9999, or caroline.jones@ivoclarvivadent.com Note: Courses scheduled are subject to change.

Join us for the Ivoclar Vivadent Welcome Reception
at THE dental event of the Year

DENTAL EXPO – INTERACTIVE

FREE to attend, ALL welcome

Register TODAY at www.dentalexpo.co.nz



2013 **dentale** PO
An Interactive Showcase

VIADUCT EVENTS CENTRE
AUCKLAND

17th - 18th May 2013



Watch your rewards grow in 2013

IvoRewards

Ivoclar Vivadent Rewards Programme



No complicated points systems, **IvoRewards** are REAL rewards for your business.

You can earn up to 15% of your annual consumable spend back in REAL rewards each year*

Call TODAY for more information on 0508 IVOCLAR

Prices do not include GST or Freight & may change without notice • Prices are not applicable to Government, Corporate or Overseas Accounts • Prices valid from 1 April 2013 – 31 May 2013
These specials are available while stocks last or until the next publication of IVOLAB • Bonus IvoRewards are only redeemable upon qualification for the IvoRewards Programme • E & OE.



AUCKLAND (Head Office)

Phone: 0-9-914 9999

Fax: 0-9-914 9990

WELLINGTON

Phone: 0-4-473 2787

Fax: 0-4-473 2656

CHRISTCHURCH

Phone: 0-3-366 0025

Fax: 0-3-377 1625

Website: www.ivoclarvivadent.co.nz

Email: orders.nz@ivoclarvivadent.com

FREephone: 0508 486 252

ivoclar
vivadent
passion vision innovation