

# Success<sup>®</sup> Injection System

## The Most Versatile Injection System in the World

Uses	Success	Ivocap <sup>®</sup>	Acron <sup>®</sup>	Palapress <sup>®</sup>	Valplast <sup>®</sup>	Flexite <sup>®</sup>
Heat-Cured	☐	■			■	
Microwave	☐		■			
Pour Resin	☐			■		
Silicone Soft Liner	☐					
Thermoplastic	☐				■	■

## Trust DENTSPLY Trubyte to develop a system you can believe in.


To develop the Success System, we consulted with the experts – you and your lab colleagues across the U.S. and Europe.

That's why you get:

- superbly fitting dentures
- precise rearticulation in the lab
- reduction in lab time on post processing adjustments
- few or no adjustments in the dental office

With the most versatile injection system in the world and premium performance from Lucitone dental resins, you can expect success every time.

For more information on making our Success part of your success, call your authorized Success Injection System dealer or call 1-800-786-0085.

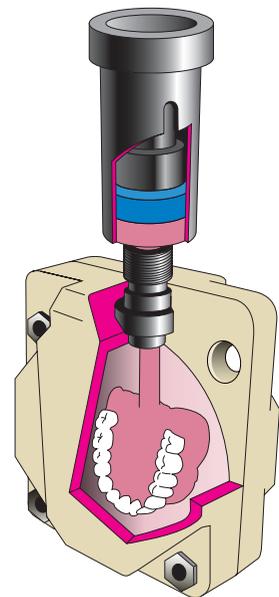

## Success<sup>®</sup> Injection System

Ivocap is a registered trademark of Ivoclar AG. ProBase is a registered trademark of Ivoclar North American, Inc. Valplast is a registered trademark of Valplast International Corporation. Flexite is a registered trademark of Rapid Injection System Corp. Success and Lucitone are registered trademarks of DENTSPLY International Inc. Listerine, Polident, Acron, Palapress, Fricke HH and St. George Formula are not trademarks of DENTSPLY International Inc.

**DENTSPLY**  
TRUBYTE

*The chemistry of a great smile.™*

DENTSPLY International, P.O. Box 872, York, PA 17405-0872 1-800-786-0085  
DENTSPLY Canada, 161 Vinyl Court, Woodbridge, ON L4L 4A3 1-800-263-1437  
www.trubyte.dentsply.com www.dentsply.com

©2002, 2004, 2006 DENTSPLY International, Inc. All right reserved.

Printed in U.S.A. Form No. 4187-B (04/06)

R<sub>x</sub> only



**DENTSPLY**  
TRUBYTE

*The chemistry of a great smile.™*

# Success<sup>®</sup> Injection System

## The Key to Success... Greater Versatility

Our high-pressure injection drives resin into a closed flask to deliver precise replication of the model. Injecting denture resin under continuous pressure fills all voids completely and compactly, eliminating the need for trial packing.

The closed-flask injection virtually eliminates the chance for change in vertical dimension of occlusion.

This kind of precision means few, if any, adjustments are necessary to properly occlude a case in the lab and deliver a great fit at chairside.

### Proven Accuracy

A study published in the *Journal of Prosthetic Dentistry*\* shows that the Success Injection System was significantly more accurate than the standard compression molding technique and saved time!


\*Journal of Prosthetic Dentistry; 1999; 82: 291-300

**Lucitone 199<sup>®</sup>** denture base and your Success Injection System can deliver lifelike aesthetics and unsurpassed impact strength. Curing is now even faster with our microwave flask.

From full dentures to flexible partials, Lucitone dental resins in your Success Injection System offer you superior performance and more versatility than ever before. Lucitone is the name to trust for high performance and true aesthetics every time.

### Lucitone 199 Makes an Impact

(Izod Impact in ft./lbs.) Data on file.

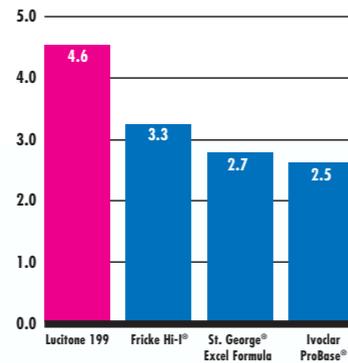

**Lucitone Fas-Por+<sup>™</sup>** pourable dental resin delivers speed to your production process without sacrificing quality. Fas-Por+ wets out faster than the leading pour acrylic. With superior color stability, Fas-Por+ creates a denture that has a smooth sheen, saving valuable finishing time.




**Luci-Sof<sup>®</sup> Injectable** denture liner is a long-term soft denture liner so durable it is a worthy addition to the Lucitone family. It's softness does not change even when soaked in harsh solutions like Listerine<sup>®</sup> and Polident<sup>®</sup>. Luci-Sof Injectable is very resistant to staining and creates an excellent bond to the denture base.



**Lucitone Clear** acrylic delivers what laboratories and dental professionals demand: strength and exceptional quality. And when you use it with the Success Injection System, you get accuracy of fit and time-saving convenience every time.

**Lucitone FRS<sup>™</sup> Flexible Dental Resin**  
Used with the Success Injection System, our flexible dental resin makes aesthetic partial dentures and very strong nightguards. The material is available in prepackaged cartridges that are heated in the Success Cartridge Furnace, then ready for precision injection.

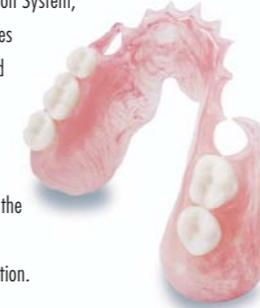

#### Benefits include:

- Easy to grind and polish
- Lightweight comfort
- Monomer and metal-free
- More impact resistant vs. conventional thermoplastics
- Lucitone 199 shades make it easy to match existing dentures
- Clear clasps – no need to match tooth shades
- UL certified furnace

### Lucitone 199 shades.



### Thermoflex<sup>™</sup> Acetal Resin

Now the Success<sup>®</sup> Injection System can inject Thermoflex acetal resin. Simply purchase a Success adapter kit. No cartridge required! Thermoflex acetal resin delivers 19 tooth colors, 3 tissue colors, excellent retention and durable performance.

More importantly, this premium-quality, pharmaceutical-grade resin brings a new dimension of simplicity and economy to your business. Now formulated in tiny pellets, Thermoflex acetal resin allows you to easily measure the precise amount of acetal resin you need. The result? Less hassle. Minimal waste.

- Aesthetic
- Flexible
- Durable
- Versatile

